



Activate Device

Endo 1

Wireless irrigation

Convenient and efficient




GUILIN WOODPECKER MEDICAL INSTRUMENT CO., LTD.



Scan and Login website
for more information

Ultrasonic irrigation – Efficiently remove residuals

With 45 kHz high frequency ultrasonic vibration, the flushing liquid can reach 3mm under the top end of irrigation tip under the action of acoustic streaming, and the chemical residue in root canal can be well cleaned.

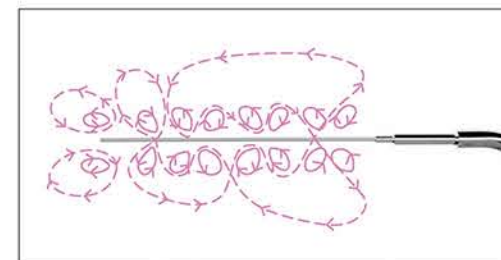

Smooth one-piece anti-break file design – Smaller amplitude and safer operation

With a small amplitude of 25 μ m, in the case of instrument fracture, the breakage is more likely to occur in the root part, which makes it easy to take out.

Upgraded irrigation file makes root canal treatment more endurable to patients.

Wireless design – Simplified operations

Eliminate the constraint of tail cord, have no need to stop operation to replace tip, and help you to complete irrigation operation with only one device.

The head of handpiece can be removed and autoclaved to avoid cross infection.

25 μ m



Handpiece can be autoclaved – Clean and reliable



Guilin Woodpecker Medical Instrument Co., Ltd.
Information Industrial Park, Guilin National High-Tech
Zone, Guilin, Guangxi, 541004 P. R. China

Tel:

Europe Sales Dept.: +86-773-5873196, +86-773-2125222

North America, South America &

Oceania Sales Dept.: +86-773-5873198, +86-773-2125123

Asia & Africa Sales Dept.: +86-773-5855350, +86-773-2125896

Fax: +86-773-5822450

E-mail: woodpecker@glwoodpecker.com, sales@glwoodpecker.com

Website: <http://www.glwoodpecker.com>

EC REP MedNet GmbH
Borkstrasse 10 · 48163 Muenster · Germany