



core white
BLEACHING LIGHT

core white

G.Comm S.r.l.

Via Donizetti, 22

20872 Cornate d'Adda (MB) Italy

T +39 039 60 60 420

F +39 039 69 26 991

info@gcomm-online.com

www.gcomm-online.com





FEATURES




TECHNIC



TECHNICAL FEATURES


Weight: approx. 800g
 Power supply: 17-24Vac
 Maximum power consumption: 30VA
 Current consumption: 1.5A
 Light pattern: 80x20mm
 Light source: 4LEDs (3.5W each)
 Emission spectrum: 420-480nm
 Maximum light power emitted: 2000mW
 Maximum distance from patient: 10cm

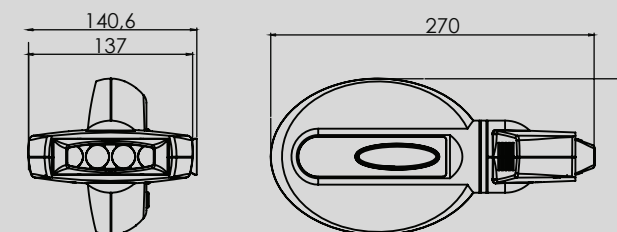

A LONG HISTORY OF QUALITY AND VERSATILITY


CoreWhite is a bleaching system for tooth whitening implemented for POLARIS or on a cart system, to grant a 360° range of motion. The cart system is robust made in steel and dust painted so as to be durable. Easy to clean does not suffer any major aggressive cleaning detergent.

LEDS SYSTEM

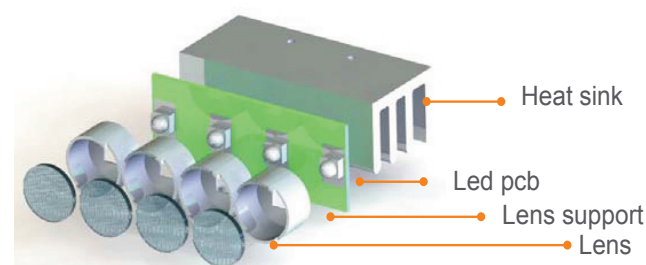


The optical system of Core White is made as shown in the figure.
 The light beams emitted by each of 4 LED sources are focused with special lenses and concentrated into one rectangular and homogeneous pattern.
 Special filters are used to eliminate I.R. on luminous pattern.

FOUR DIFFERENT WORKING PROGRAMS

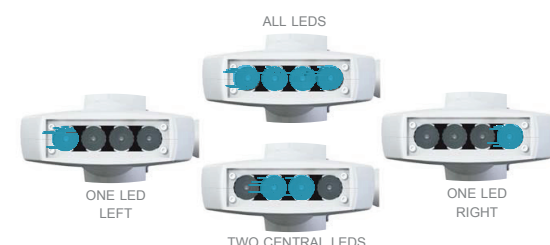

1-20min TREATMENT TIME ADJUSTABLE

Cart application

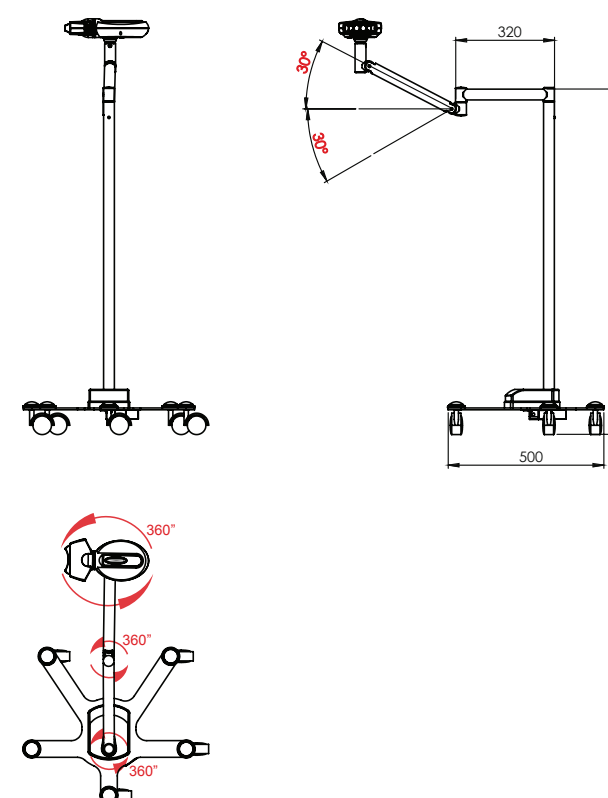

Polaris application

